


TGIRT ssDNA-seq


Douglas Wu

Overview

Procedures

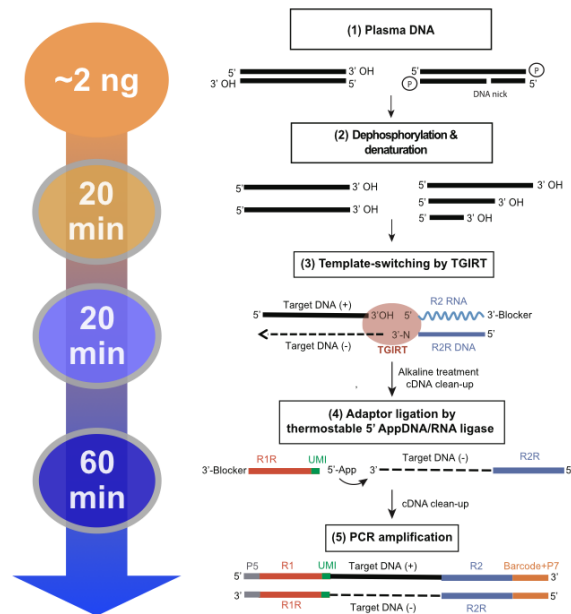

Final sequence structure

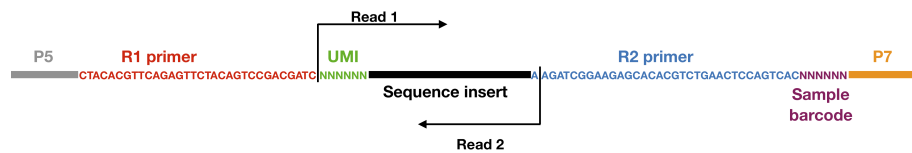

Reagents for TGIRT-DNA-seq

Buffers

You can scale up accordingly.

DNA TS Buffer I

| Material | Volumn | Final concentration |
|----------------------|----------|---------------------|
| 2M MgCl ₂ | 105 uL | 21 mM |
| 2M Tris-HCl pH 7.5 | 430 uL | 86 mM |
| H ₂ O | 9.465 mL | |
| Total | 10 mL | |

DNA TS Buffer II (420 mM)

| Material | Volumn | Final concentration |
|----------------------|-----------|---------------------|
| 2.4M NaCl | 8.4 mL | 200uM |
| 2M MgCl ₂ | 87.5 uL | 17.5 mM |
| 2M Tris-HCl pH 7.5 | 350 uL | 70 mM |
| H ₂ O | 1162.5 uL | |
| Total | 10 mL | |

TE Buffer

| Material | Volumn | Final concentration |
|---------------------|---------|---------------------|
| 2 M Tris HCl pH 7.5 | 50 uL | 10 mM |
| 0.5 M EDTA pH 8 | 20 uL | 1 mM |
| H ₂ O | 9.93 mL | |
| Total | 10 mL | |

Primers

DNA primers

- R2R (TTN): GTG ACT GGA GTT CAG ACG TGT GCT CTT CCG
ATC TTN
- 6N UMI R1R: /5'Phos/ NNN NNN GAT CGT CGG ACT GTA GAA
CTC TGA ACG TGT AG /3'C3Sp/
- Illumina Barcoded PCR primer: CAA GCA GAA GAC GGC ATA CGA
GAT NNNNNN GTG ACT GGA GTT CAG ACG TGT GCT CTT CCG
ATC T
- Multiplex PCR primer: AAT GAT ACG GCG ACC ACC GAG AT C
TAC ACG TTC AGA GTT CTA CAG TCC GAC GAT C

RNA primer

- R2 (TTN): AAG AUC GGA AGA GCA CAC GUC UGA ACU CCA GUC
AC/3'C3Sp/

Preparation of 5'App-R1R

| Material | Volumn |
|------------------|--------|
| H ₂ O | 13 uL |
| 10X rxn buffer | 2 uL |
| 1mM ATP | 2 uL |
| 10 uM 6N UMI R1R | 1 uL |
| Mth RNA ligase | 2 uL |
| Total | 20 uL |

- Incubate at 65°C for 1 hr
- Incubate at 85°C for 5 min
- Cleanup using NucleoSpin Gel and PCR Clean-up (use Binding buffer
NTC, Macherey-NagelTM 740654.100)

Protocol

Pre-anneal template-switching duplex

| Material | Volumn |
|-----------------|--------|
| TE buffer | 12 uL |
| 10uM R2R primer | 1.5 uL |
| 10uM R2 primer | 1.5 uL |
| Total | 15 uL |

* 15 uL is enough for 7 rxns, can scale up if needed

- Incubate at 82°C for 2 min
- Cool down to 25°C with 10% ramp/0.1°C/second

Prepare DNA template

| Material | Volumn |
|-----------------|---------|
| DNA TS Buffer I | 3 uL |
| 0.5M DTT | 0.2 uL |
| DNA sample | 9 uL |
| FastAP | 1 uL |
| Total | 13.2 uL |

- Incubate at 37°C for 20 min (start preparing template-switch complex)
- Denature at 95°C for 3 min

Prepare template-switch complex

| Material | Volumn |
|---------------------|--------|
| DNA TS Buffer II | 2 uL |
| Pre-annealed duplex | 2 uL |
| TGIRT | 2 uL |
| Total | 6 uL |

- Incubate at room temperature for 30 min

R2-adapter addition and cDNA synthesis

- Add 6 uL template-switch complex into 13.2 uL of DNA template
- Add 0.8 uL dNTP (25 mM) to the mixture
- Incubate at 60°C for 20 min
- Add 1 ul of 5 M NaOH and incubate at 95°C for 3 min
- Add 1 ul 5 M HCl to neutralize pH
- Clean-up with Qiagen MinElute Reaction Cleanup Kit

R1-adapter ligation

| Material | Volume |
|---------------------------|--------|
| 10X NEB-1 buffer | 2 uL |
| 50mM MnCl ₂ | 2 uL |
| 10uM 5'App-R1R | 4 uL |
| cDNA | 10 uL |
| 5'App thermostable ligase | 2 uL |
| Total | 20 uL |

- Incubate at 65°C for 1 hr
- Incubate at 90°C for 3 min
- Clean-up with Qiagen MinElute Reaction Cleanup Kit

PCR amplification and addition of sequencing primers

| Material | Volume |
|--------------------------|--------|
| Barcoded PCR primer | 2.5 uL |
| Multiplex PCR primer | 2.5 uL |
| Ligation product | 20 uL |
| KAPA hotstart master-mix | 25 uL |
| Total | 50 uL |

- PCR cycles:
 - 98°C 30 sec, 1 cycle
 - 98°C 45 sec, 60°C 15 sec, 72°C 30 sec 11 cycles
 - 72°C for 5 min, 1 cycle
 - hold at 4°C

Final cleanup

- Add 65 uL Agencourt AMPure XP beads into the PCR product and pipette the mixture to a 2 mL microcentrifuge tube
- Incubate at room temperature for 10 min
- Put the tube on magnetic rack and wait 5 min
- Remove all liquid
- Do twice: Add 200 uL 80% EtOH, wait 30 sec and remove
- Quick spin down, and put the tube on magnetic rack, remove residual EtOH using 20P pipette
- Add 31 uL H₂O, mix well and incubate in room temperature for 10 min
- Put the tube on magnetic rack and wait 5 min
- Pipette 30 uL of elute to a clean tube